

IMPROVING THE LIVES OF COUNTLESS CHILDREN AND ADULTS

ABT-Advanced Brain Technologies, the creators of iListen, are the developers of innovative brain training products that maximize human potential.

Our products are improving the lives of countless children and adults. By training brain function at its core, it allows us to address diverse issues such as memory, listening, sensory processing, brain longevity, peak-performance and more.

We have organized an expanding global team of multidisciplinary specialists applying over 30 years of clinical experience. The ABT team combines advances in neuroscience, psychoacoustic and music research with recent developments in the field of brain plasticity, the brain's natural ability to change itself.

ABT created its exclusive recording label, ABT Music®, to focus on the company's unique approach toward researching, selecting, arranging, recording, editing, mastering and clinically testing high-quality acoustic music. ABT Music invests continually in research and development at its state-of-the-art production facilities to produce award-winning music for the benefit and enjoyment of listeners.

Current products include professionally administered therapeutic products: The Listening Program®, Spatial Surround® HD, and iListen™, as well as interactive software and music programs: BrainBuilder®, Sound Health®, Music for Babies™, and Shelly's Interactive Computer Learning Series™.

 **ABT** ADVANCED BRAIN
TECHNOLOGIES
5748 South Adams Avenue Parkway
Ogden, Utah 84405

iListen™

iListen REGISTRATION

To receive service or warranty support, you must register your iListen. Go to thelisteningprogram.com and click on the "TLP Listener Resources" button located on the home page. Then select "iListen Registration" and complete the form.

WARRANTY

ABT-Advanced Brain Technologies warrants the bone conduction amplifier and bone conduction headphones (Product) to be free from defects in material and workmanship for a period of 90 days. If a Product proves to be defective in material or workmanship during the warranty period, ABT will, at its sole option, repair or replace the Product with a similar Product. Replacement Product or parts may include remanufactured or refurbished parts or components. The replacement unit will be covered by the balance of the time remaining on the customer's original limited warranty. Replacement of lost files due to illegal file transfer attempts is specifically excluded under this warranty. The iPod® includes a one year warranty directly from Apple, Inc.

DISCLAIMER

The programs on iListen have not been reviewed by the Food & Drug Administration (FDA). The programs on iListen are not intended to diagnose, treat, cure, or prevent any disease.

©2007 ABT-Advanced Brain Technologies | All Rights Reserved
iListen™ is a trademark owned by Advanced Brain Technologies. iPod® is a registered trademark of Apple, Inc. Packaging, system and product configuration, ©2007 Advanced Brain Technologies, Ogden, Utah 84405. All rights reserved.

FOR ADDITIONAL INFORMATION OR SUPPORT, PLEASE CONTACT YOUR PROVIDER:

Visit us online at www.thelisteningprogram.com

 **ABT** ADVANCED BRAIN
TECHNOLOGIES

iListen™
Quick Start Guide

CONNECTING YOUR UNIT



- 1 To connect iPod to Headphone Amplifier: Connect the enclosed 1 ft. cable from the iPod Headphones port to the port marked IN on the Headphone Amplifier.
- 2 To connect Headphones to Headphone Amplifier: Connect the ALL BLACK connector on the Headphone cable to the port marked BC on the Headphone Amplifier. Then connect the connector with the BLUE RING on the Headphone cable to the port marked AC on the Headphone Amplifier.
- 3 To power Headphone Amplifier ON/OFF: Toggle switch UP to power on and DOWN to power off. PWR indicator light will glow blue when unit is on.

PROTECT HEARING! To protect your hearing health, the maximum volume level of the iPod has been pre-set and locked by ABT. To change this setting, please contact your Provider.

www.thelisteningprogram.com


PLAYING A PROGRAM

1. Unit is connected. DO NOT place Headphones on at this time.
2. Push MENU on the iPod Click Wheel. Make sure iPod nano is not on HOLD.
3. Select MUSIC. Choose your selection by lightly moving your thumb around the Click Wheel and clicking the Center Button on the iPod.
4. Select ALBUMS and click the Center Button on the iPod.
5. Select a PROGRAM and click the Center Button on the iPod.
6. Adjust iPod volume to HIGH.
7. Now power on Headphone Amplifier. On back of the Headphone Amplifier, place the BC volume setting, AC L and AC R levels all at "0".
8. Place Headphones on. Select a MODULE and click Center Button on the iPod. Module will begin to play.
9. On the back of the Headphone Amplifier, slowly adjust the BC volume control and the AC L and AC R levels given to you by your Provider.
10. Be sure the volume is at a comfortable level. For any questions about volume settings, please contact your Provider.

WARNING! DO NOT attempt to connect the iPod to your computer or file loss may result. Neither Advanced Brain Technologies nor its Providers is responsible for loss of files. This unit is equipped with copy protection. Do not attempt to copy files from iListen.

NOTE: To re-charge the iPod, please use the AC Power Recharger included in the iListen package. The Headphone Amplifier uses two (2) AA batteries.

iPod nano

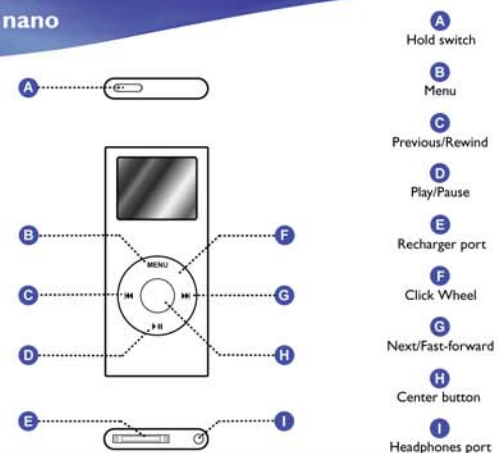

TO	DO THIS
Turn on iPod nano	Press any button.
Turn off iPod nano	Press and hold Play/Pause.
Turn on the backlight	Press any button or use the Click Wheel.
Go directly to the main menu	Press and hold Menu.
Change the volume	From the Now Playing screen, move your thumb around the Click Wheel.
Place iPod nano on HOLD (Disable iPod buttons to prevent accidental volume increases and control function mishaps while listening)	Slide the Hold switch to HOLD. An orange strip will appear when the iPod is on HOLD.
Re-charge iPod nano	Connect AC Power Recharger.
Reset, "unfreeze" an iPod nano that does not respond, won't turn on, or remains blank	Check the Hold switch to make sure it's not on HOLD. Check your battery to make sure it is charged. Slide the Hold switch to HOLD and back again. Then press the Menu and Center buttons at the same time until the Apple logo appears, about 6-10 seconds.
Service your iPod nano	Complete the online iPod request form found by selecting the Support tab at the top of apple.com . You can also take your iPod to an Apple retail store or an Apple Authorized Service Provider.